



NA4 Glycan

NA4 N-linked oligosaccharide

Catalog Number

SA-NA4

Description

Asialo-, galactosylated, tetra-antennary complex-type N-glycan (oligosaccharide).

Sources

Found on several mammalian glycoproteins and human alpha-1 acid glycoprotein (hAAGP). This product is typically purified from the oligosaccharide pool released from hAAGP by hydrazinolysis using a combination of HPLC and glycosidase digestion.

Form

Dry. Dried by centrifugal evaporation from an aqueous solution.

Molecular Weight


2373

Purity

> 90% pure as assessed by a combination of 1H-NMR and HPLC.

Shipping

The product can be shipped at ambient when dry. After dissolution, ship on dry ice.

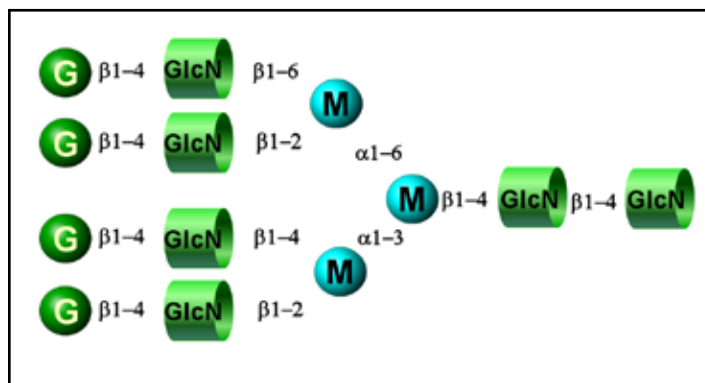

Storage

Keep at -20°C both before and after dissolution. This product is stable for at least 5 years as supplied.

Handling

Allow the unopened vial to reach ambient temperature and tap unopened on a solid surface to ensure that most of the lyophilized material is at the bottom of the vial. Gently remove the cap, add the desired volume of reconstitution medium, re-cap and mix thoroughly to bring all the oligosaccharide into solution. For maximal recovery of oligosaccharide, ensure that the cap lining is also rinsed and centrifuge the reconstituted vial briefly before use. Ensure that any glass, plasticware or solvents used are free of glycosidases and environmental carbohydrates.

Safety

This product is non-hazardous and has been purified from natural sources certified to be free of all hazardous material including pathogenic biological agents.

Certificate of Analysis

SA-NA4 NA4 Glycan
Specifications


Product NA4 Glycan

Catalog no. SA-NA4

Lot no. A675-03

Purity > 90% pure as assessed by a combination of $^1\text{H-NMR}$ (see Fig 1) and HPLC.

Figure 1 : 500 MHz $^1\text{H-NMR}$ of NA4 Glycan (Cat. No.SA-NA4, Lot No. A675-03t)

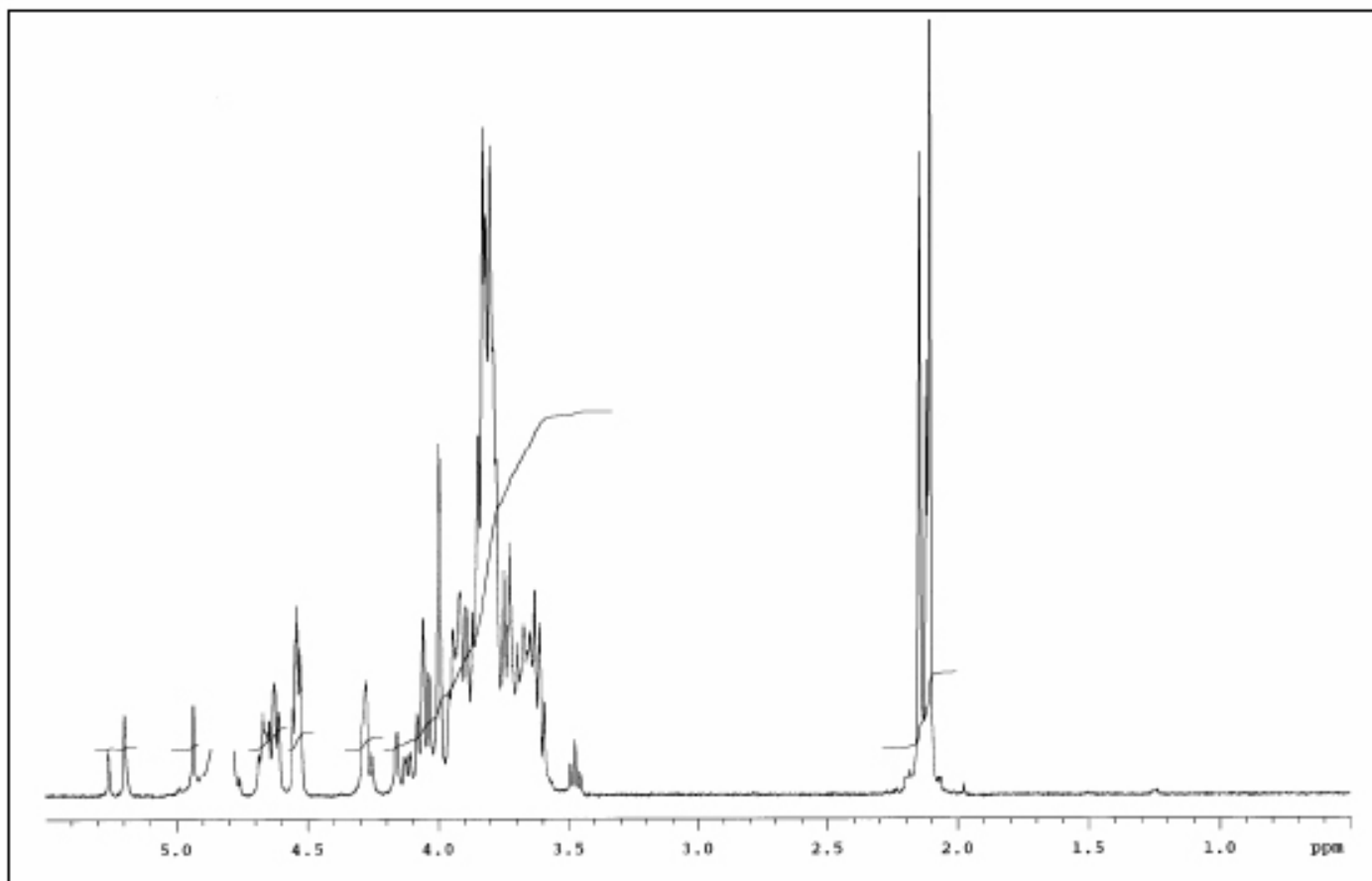

Warranties and liabilities

QA-Bio, LLC warrants that the above product conforms to the specifications described herein. Should the product fail for reasons other than through misuse QA-Bio, LLC will, at its option, replace free of charge or refund the purchase price. This warranty is exclusive and QA-Bio, LLC makes no other warrants, expressed or implied, including any implied conditions or warranties of merchantability or fitness for any particular purpose. QA-Bio, LLC shall not be liable for any incidental, consequential or contingent damages.

This product is intended for *in vitro* research only.

updated 9/23/201